

POWER Q-6000 PLUS

Professional device for compressive limbs therapy with 6 compression chambers and 2 outputs. With POWER Q-6000 Plus is possible to customize each treatment according to needs.

POWER Q-6000 Plus promotes limb blood flow by continuously repeating swelling and contraction. Moreover it promotes the excess interstitial fluids elimination and thus aiding the correct function of the lymphatic system.



TECHNICAL FEATURES

- ✓ Power supply 220-240 V 50/60 Hz
- ✓ Power consumption 80 W
- ✓ Adjustable time 0÷99 minutes
- ✓ Adjustable pressure 0÷240 mmHg ±10%
- ✓ 2 output channels (up to 3 connectable applicators)
- ✓ 5 programs
- ✓ 8 air speed levels
- ✓ LCD Graphic display
- ✓ Weight 9.5 Kg

MEDICAL DEVICE CLASS IIa

POWER Q-6000 Plus is a medical device CE0476

MAIN APPLICATIONS

Lymphoedema treatment

Limb paralysis

Passive compression for venous disorders

Rheumatoid arthritis


Arm drainage after mastectomy

Treating hematomas after surgeries

Anti-cellulite treatment in the initial stages

Post operation prophylaxis in venous thrombosis surgery

PROGRAMS

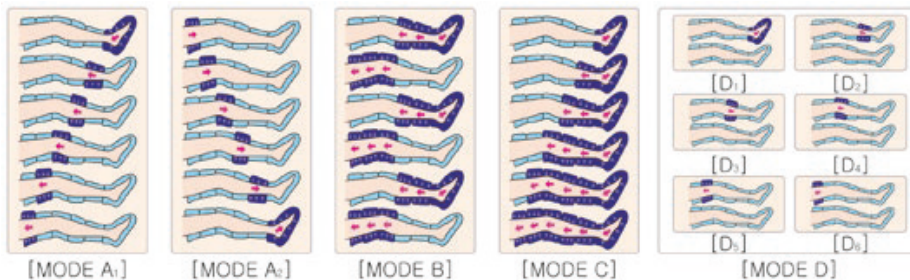

DIMENSIONS

420 x 230 x 310 mm

POWER Q-6000 PLUS

HOW DOES IT WORK

Compressive limbs therapy is able to promote a right venous circulation by decreasing muscle tension caused by stress or chronic and acute pains: thanks to its pump action it favors venous return by increasing tissues blood circulation and finally their physiological renewal.

Compressive limbs therapy can solve blood stagnation problem in damaged sanguine vessels or in body zones not correctly bedewed.

Most common applications are:

- ✓ Body massage increase
- ✓ Blood flow promotion thanks to a slow and progressive pressure
- ✓ Lymphatic drainage
- ✓ Lower limbs relief and comfort
- ✓ Substances stagnation prevention in peripheral tissues
- ✓ Blood insufficiency prevention caused by an inadequate pump action
- ✓ Muscle tone loss in disabled or paralyzed patients. Compressive limbs therapy substitutes the manual massage for trophic purposes.

POWER Q-6000 Plus is a professional device for compressive limbs treatments. It promotes limb blood flow by continuously repeating swelling and contraction. The external compression promotes the excess interstitial fluids elimination and thus aiding the correct function of the lymphatic system.

POWER Q-6000 Plus offers 5 different inflation modes of the air chambers, 8 speed levels of the air flow inside the chambers and the ability to set the time and the pressure of the therapy. In this way it is possible to customize each treatment according to the therapeutic needs of each patient.

